



SAFElife™

Multi-Drug Urine Test Cup

- Easy to use
- Customizable panels
- Rapid & Convenient
- Reliable Results

The most comprehensive panel of **24** CLIA waived drug tests.

Now available:

6-MAM 10 ng/mL, Fentanyl 1 ng/mL, Norfentanyl 5 ng/mL, THC 20 ng/mL and Tramadol 100 ng/mL.



6-Monoacetylmorphine (6-MAM): 10 ng/mL
Amphetamines (AMP): 500 or 1000 ng/mL
Secobarbital (BAR): 300 ng/mL
Buprenorphine (BUP): 10 ng/mL
Oxazepam (BZO): 300 ng/mL
Cocaine (COC): 150 or 300 ng/mL
EDDP: 300 ng/mL

Fentanyl (FTY): 1 ng/mL
Methylenedioxymethamphetamine (MDMA): 500 ng/mL
Methamphetamine (MET/mAMP): 500 or 1000 ng/mL
Morphine (MOP): 300 ng/mL
Methadone (MTD): 300 ng/mL
Opiates (OPI): 2000 ng/mL
Norfentanyl (NFTY): 5 ng/mL

Oxycodone (OXY): 100 ng/mL
Phencyclidine (PCP): 25 ng/mL
Propoxyphene (PPX): 300 ng/mL
Nortriptyline (TCA): 1000 ng/mL
Cannabinoids (THC): 20 or 50 ng/mL
Tramadol (TRA): 100 ng/mL

Quick Test Process

Please refer to the Package Insert for detailed process information.


- 1 Remove the test cup lid and deposit urine directly into the test cup until it reaches the Minimum Fill Level.



- 2 Peel off the label from right to left to view the results.



- 3 Read the adulteration test results immediately, at 30 seconds and at 45 seconds. Read the drug test results at 5 minutes.

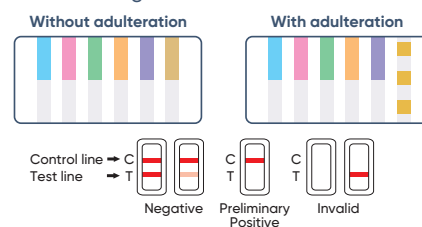

Ordering Information

Item Code	Description	Classification
CDOA-L11135	13-Panel T-Cup Compact; AMP500,BAR300,BUP10,BZO300,COC150,FTY1,MDMA500,MET500,MOP300,MTD300,OXY100,PCP25,THC50	CLIA waived
CDOA-L11145	14-Panel T-Cup Compact; AMP500,BAR300,BUP10,BZO300,COC150,FTY1,MDMA500,MET500,MOP300,MTD300,OXY100,PCP25,TCA1000,THC50	CLIA waived
CDOA-L11145A3	14-Panel T-Cup Compact; AMP500,BAR300,BUP10,BZO300,COC150,FTY1,MDMA500,MET500,MOP300,MTD300,OXY100,PCP25,TCA1000,THC50 with OX/pH/S.G.	CLIA waived
CDOA-L12145A3	14-Panel T-Cup Compact; AMP500,BAR300,BUP10,BZO300,COC150,FTY1,MDMA500,MET500,MOP300,MTD300,OXY100,PCP25,THC50,TRA100 with OX/pH/S.G.	CLIA waived

Rel:
2025/12/04



+1 (630) 468-2199 || wondfo@wondfousa.com || WondfoUSA.com

Wondfo® || bringing science to life